

# TANZANITE (AMC 2520)



An electrically powered 60° tilt table suitable for various applications and customisable to suit your needs.

Uniquely designed to raise in height when in the tilt position to allow practitioners easier access to patients.

Designed to reduce OH&S risks and minimise manual handling.



**5 Year warranty**  
on electrical components



**Lifetime warranty**  
on frame and workmanship



**Australian Made**

## Features

Adjustments of height (hi/lo), back rest and tilt made with high quality electric actuators, reducing OH&S risks.

Designed with a manually operated foot rest to support patients during tilt.

Supplied with hand control switch.

Standard with battery pack and central wheel locking.

Low minimum height gives patients the ability to easily get on and off the couch.

Carefully considered base design allows easy side access for practitioners and patients.

Access under couch for patient lifting devices.

125mm lockable wheels allow easy movement of the couch.

## Basic specifications

Width: 610mm

Length: 2000mm

Couch tilt: 60°

Back rest tilt: 70°

Height adjustment: 550mm -1000mm

Lifting capacity: 200kg

Wheels: 125mm central locking

## Optional customisations

Split foot rest

Cardio drop section with gas strut

Cardio cut out on left or right side

Foot switch control for hands free operation

Head cushion

Paper roll holder

Wide variety of fabric colours available


Padded swivel arm rest set

Standard arm rest set


Fold down side rails

# TANZANITE (AMC 2520)

## Side view

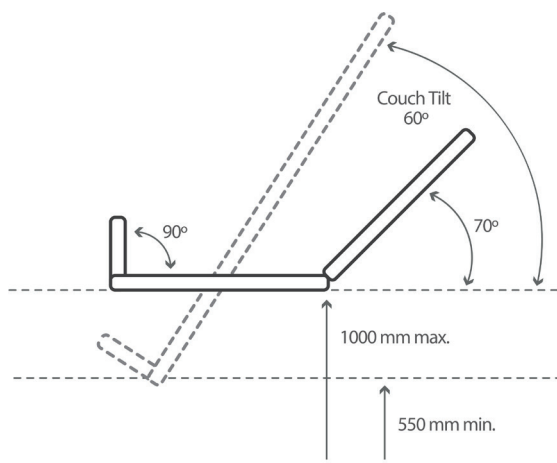

## Top view

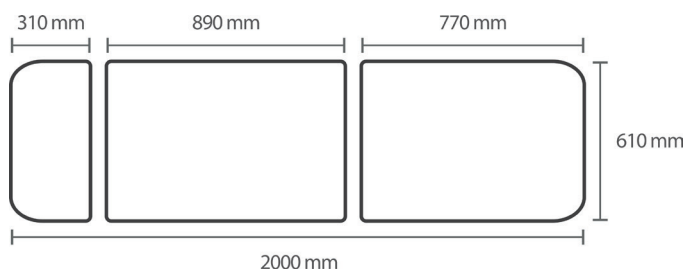

## Optional customisations



Cardio drop section



Cardio cut out



Fold down side rails



Paper roll holder

## Detailed specifications

Head section length: 770mm

Middle section length: 890mm

Foot section length: 310mm

Vinyl specifications:

UV resistant

Flame resistant

Easy to clean

Reduces growth of bacteria

Complies with AS 1530.3 standards

Variety of colours available

Foam specifications: 50mm foam

Frame material and colour:

Powder coated white frame

Australian steel

Wheel system:

Central locking castors

Drive actuators:

Couch tilt: Linak electric actuator system

Back rest: Linak electric actuator system

Hi/lo: Linak electric actuator system

Electrical specifications: 230V 50-60 Hz

Electrical standards: IPX4 rating

Lift capacity: 200kg

Shipping weight: 70 - 85 kg

Shipping dimensions:

2130mm x 820mm x 700mm

Freight charges:

Not included in sale price.

Contact our sales representatives for a quote.

ARTG number: 208309

Warranty:

Lifetime on frame and workmanship

5 years on electrical components

Suitable for use in:	GENERAL EXAMINATION	ULTRASOUND	TREATMENT ROOMS	GYNAECOLOGY & OBSTETRICS	VASCULAR	CARDIOLOGY	DERMATOLOGY	ENT	PODIATRY	COSMETIC SURGERY	RENAL DIALYSIS	CHEMO / DAY ONCOLOGY	BLOOD COLLECTION	DAY RECOVERY	ICU	DAY SURGERY	SPECIALIST APPLICATIONS