

TOPAZ [AMC 2120]



Three section treatment couch that allows you to get closer to your patients.

Electric adjustments of height, back rest and each individual leg rest.

Designed to reduce OH&S risks and minimise manual handling.



5 Year warranty
on electrical components



Lifetime warranty
on frame and workmanship



Australian Made

Features

Movements powered by high quality electric actuators, reducing OH&S risks.

Supplied with hand control switch and two arm rests.

Independently moving leg rests.

Low minimum height gives patients the ability to easily get on and off the couch, including wheelchair transfer.

Carefully considered base design allows easy side access for patients and practitioners.

Access under couch for patient lifting devices.

125mm lockable wheels allow effortless movement of the couch.

Basic specifications

Width: 610mm

Length: 1870mm

Back rest tilt: 75°

Foot rest tilt: 45°

Height adjustment: 550mm -1000mm

Lifting capacity: 200kg

Wheels: 125mm lockable castors

Optional customisations

Extra width - 660mm

Foot switch control, for hands free operation

White LED light, attached to the couch or to a moveable stand

Battery pack

Head cushion

Paper roll holder

Wide variety of fabric colours available


Padded swivel arm rest set, suitable for Phlebotomy

Back cushion


Heine light

TOPAZ [AMC 2120]

Side view

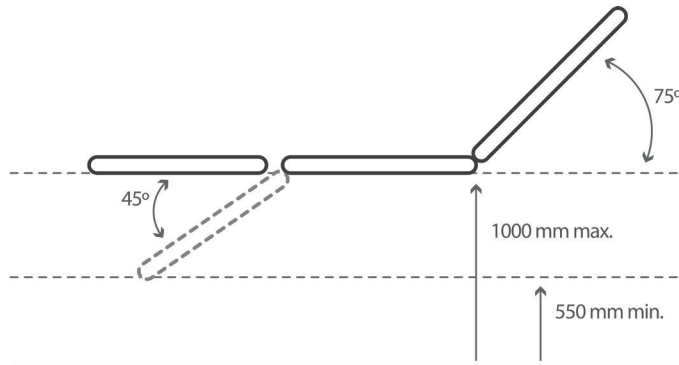

Top view

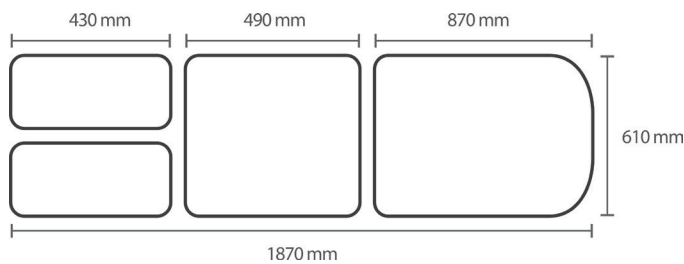

Optional customisations



White LED light



Head cushion



Padded, swivel arm rest



Paper roll holder

Detailed specifications

Head section length: 870mm

Middle section length: 490mm

Foot section length: 430mm

Vinyl specifications:

UV resistant

Flame resistant

Easy to clean

Sanitised to inhibit growth of bacteria

Complies with AS 1530.3 standards

Variety of colours available

Foam specifications: 50mm foam

Frame material and colour:

Powder coated white frame

Australian steel

Wheel system:

Individual 125mm brake castors made by Tente

Drive actuators:

Hi/lo: Dewert electric actuator system

Back rest: Dewert electric actuator system

Right foot rest: Dewert electric actuator system

Left foot rest: Dewert electric actuator system

Electrical specifications: 240V 50H 26DVC 105W

Electrical standards: IPX4 rating

Lift capacity: 200kg

Shipping weight: 70 - 85 kg

Shipping dimensions:

2130mm x 820mm x 700mm

Freight charges:

Not included in sale price.

Contact our sales representatives for a quote.

ARTG number: 208309

Warranty:

Lifetime on frame and workmanship

5 years on electrical components

| | | |
|----------------------|--|---|
| Suitable for use in: | <ul style="list-style-type: none"> ● GENERAL EXAMINATION ● ULTRASOUND ● TREATMENT ROOMS ● GYNAECOLOGY & OBSTETRICS ● VASCULAR ● CARDIOLOGY ● DERMATOLOGY ● ENT ● PODIATRY ● COSMETIC SURGERY | <ul style="list-style-type: none"> ● RENAL DIALYSIS ● CHEMO / DAY ONCOLOGY ● BLOOD COLLECTION ● DAY RECOVERY ● ICU ● DAY SURGERY ● SPECIALIST APPLICATIONS |
| | | |